


# Cavafix® Duo


Double-lumen venous catheter utilizing the Braunula® principle

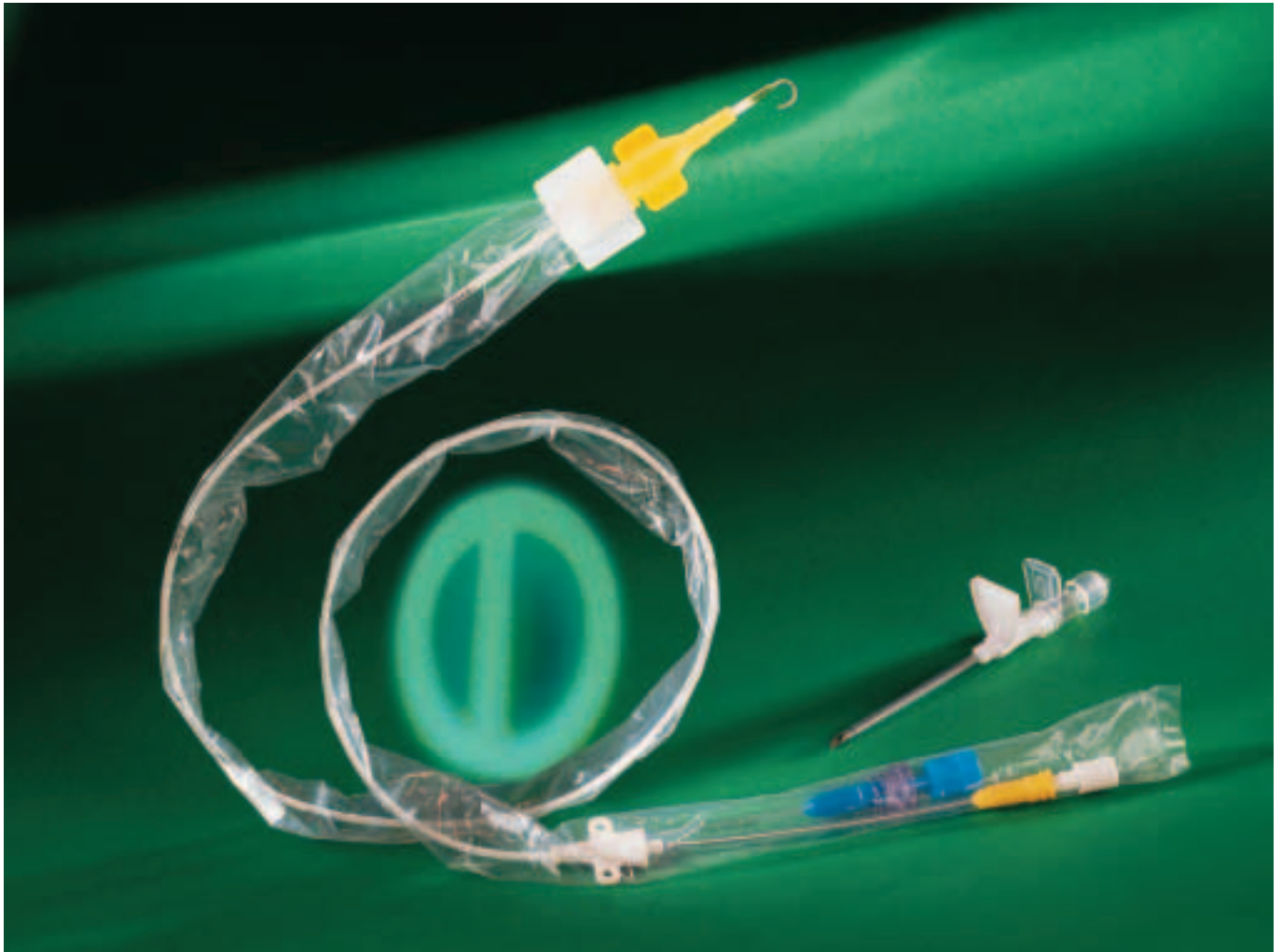

- Avoidance of incompatibilities with two separate channels for drug application
- Reduction of puncture-related infections with a single vessel access
- Rapid and reliable catheter placement under simultaneous ECG lead

# Cavafix® Duo

Double-lumen venous catheter utilizing the Braunula® principle

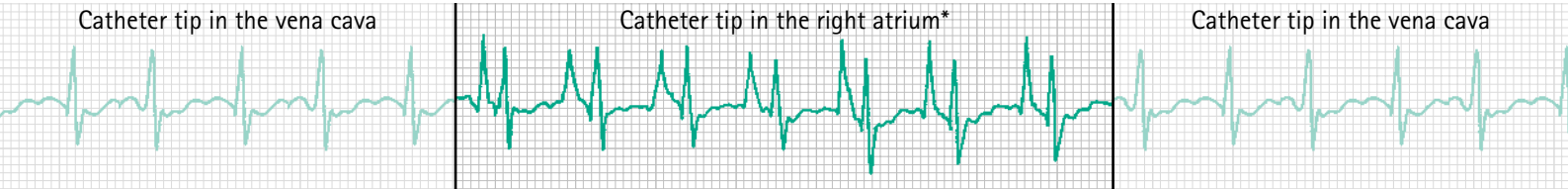

## Exact catheter positioning without x-ray

- Simultaneous monitoring of catheter position via right-atrium ECG lead allows for immediate correction of catheter position in the event of incorrect placement
- P-wave changes in the ECG curve indicate the proper catheter position\*




## Quick connection to ECG monitor

- The ECG connection cable guarantees safe and secure impulse transmission to the monitor
- The Certodyn® universal adapter can be combined with nearly all patient cables
- Switching between external and internal ECG lead is simple and safe

## Safsite® – The safety connector

- Opens and closes during connection or disconnection of the luer attachment
- Reliably prevents blood or other liquids from escaping, air is kept out of the system

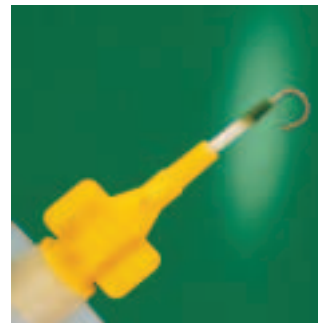

## Catheter with Soft-Tip and guidance stylet

- The Soft-Tip reduces vein perforations and injuries of the intima
- The guide stylet with its highly flexible J-tip facilitates problem-free advancement of the catheter

## Presentations

Product name	Code	Ø mm	Catheter G	Length cm	
Cavafix® Duo 435	415 1801	2,1	16/18	32	Vena jugularis
Cavafix® Duo 455	415 1810	2,1	16/18	45	Vena jugularis
Cavafix® Duo 458	415 1828	2,1	16/18	45	Vena subclavia
Cavafix® Duo 475	415 1836	2,1	16/18	70	Vena basilica
Cavafix® Duo 225	415 1453	1,4	18/20	20	Vena jugularis
Cavafix® Duo 255	415 1461	1,4	18/20	45	Vena jugularis
Cavafix® Duo 257	415 1470	1,4	18/20	45	Vena subclavia
Cavafix® Duo 265	415 1488	1,4	18/20	60	Vena basilica
Certodyn® universal adapter	415 0228				
Certodyn® universal adapter Paed	415 0724				

**B | BRAUN**  
HOSPITAL CARE

B. Braun Melsungen AG  
P. O. Box 11 20  
D-34209 Melsungen  
Tel +49 (0) 56 61 71-0  
Fax +49 (0) 56 61 71-45 55  
www.bbraun.com