
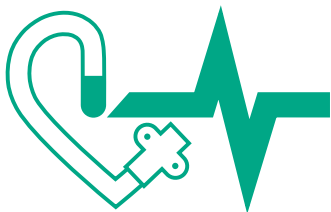


# Certofix<sup>®</sup> Quinto

For more safety  
with poly-traumatised patients

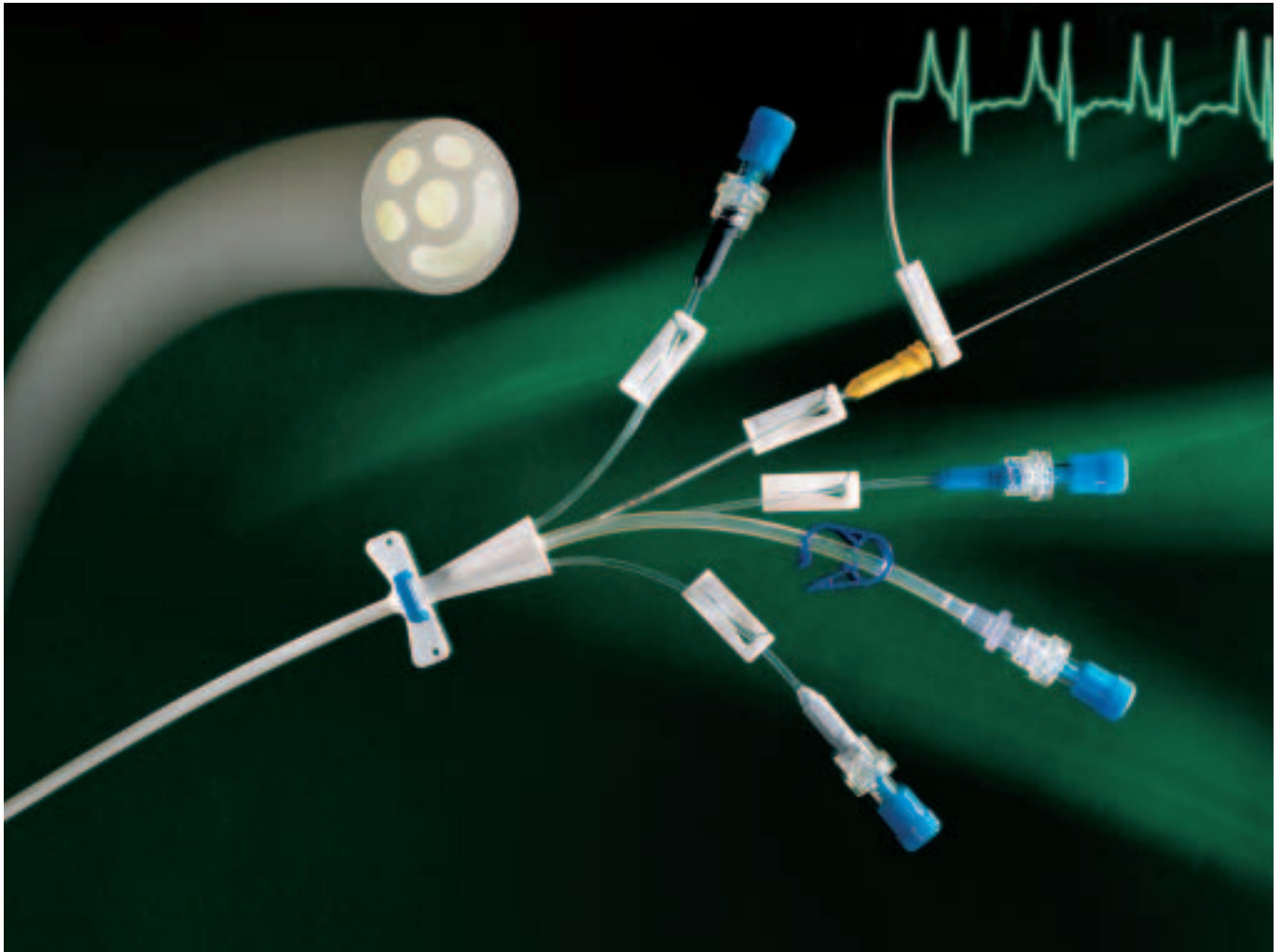

Five-lumen central venous catheters  
for placement using Seldinger technique  
with simultaneous ECG lead option

# Certofix<sup>®</sup> Quinto

For more safety  
with poly-traumatised patients

## Five-lumen catheter


- for poly-traumatised patients and for vascular and abdominal surgery
- for an optimal intensive care with complex infusion therapies



- for blood transfusions in case of extensive volume deficiency
- for the avoidance of drug incompatibilities

## Catheter with soft tip

- reduces vein perforations and intima lesions with the flexible, soft polyurethane tip

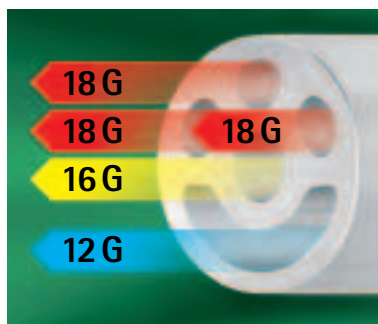


## Optimized section

- 18 G channels for the parallel application of e.g. catecholamines, barbiturates, diuretics etc.
- 16 G channel for the parenteral alimentation
- 12 G channel for transfusions or volume fluid substitution

## Simultaneous positioning check during catheterization via the atrial ECG

- enables immediate corrective measures to be made in case of initial incorrect positioning of the catheter
- helps prevent subsequent complications, particularly the risk of pericardiac perforation and cardiac arrhythmias

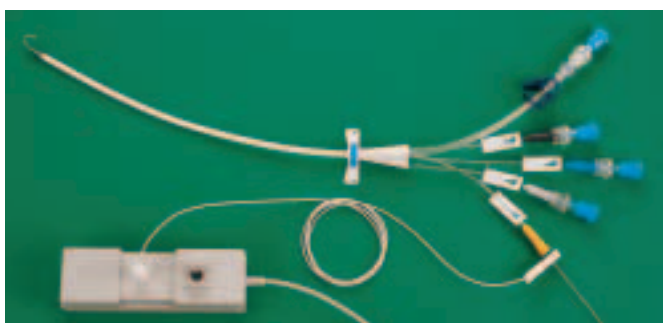

- no X-ray required, no radiation exposure
- no incidents with contrast medium
- effortless identification of the P wave and thus clear verification of position
- no subsequent recatheterization due to incorrect positioning
- reduced costs of equipment and staff

## Product Identification

Name	Code number	Catheter/Lumen Ø G	Length cm
Certofix <sup>®</sup> Quinto SB 1215	0416 6833	16/18/18/18/12	15 cm
Certofix <sup>®</sup> Quinto V 1215	0416 6841	16/18/18/18/12	15 cm
Certofix <sup>®</sup> Quinto SB 1220	0416 6850	16/18/18/18/12	20 cm
Certofix <sup>®</sup> Quinto V 1220	0416 6868	16/18/18/18/12	20 cm
Certofix <sup>®</sup> Quinto SB 1230	0416 6876	16/18/18/18/12	30 cm
Certofix <sup>®</sup> Quinto V 1230	0416 6884	16/18/18/18/12	30 cm

**B | BRAUN**  
HOSPITAL CARE

**B. Braun Melsungen AG**  
P. O. Box 11 20  
D-34209 Melsungen  
Tel (0 56 61) 71-0  
Fax (0 56 61) 71-45 55  
www.bbraun.com